

MDF-86V340E

-86°C Ultra Low Temperature Freezer




340 liters
19,200 samples



PROMED AUTO-CASCADE SYSTEM SERIES

-86°C ULT freezer is specially designed and manufactured for long term storage of various biological products, including viruses, germs, erythrocytes, leucocyte and cutis. Applications can be found in blood banks, hospitals, edipemic prevention service, and research institutes, laboratories in electronic and chemical plants and biological engineering institutes.

PROMED AUTO-CASCADE SYSTEM

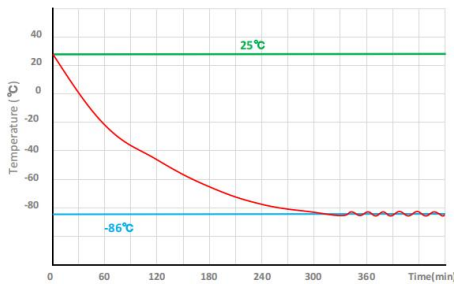

Pull down test at ambient temp.




Advanced thermal insulation

The insulation are designed with the latest generation of VIP technology to achieve thermal conductivity of 0.005 W(m.K). That makes ULTs more efficient in Insulation and less heat loss, so as to achieve the purpose of energy saving as well as space-saving.



Super-cooling technology

The powerful cooling system drives the unique mixed refrigeration system for fast cooling down. This system is stable, energy saving and easy maintenance.

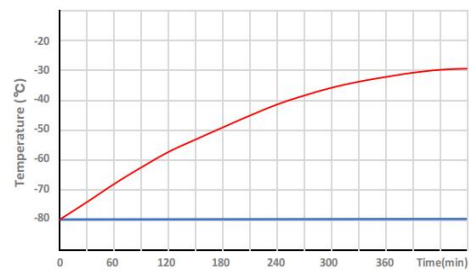

Pull up test at ambient temp. 25°C



MDF-86V588E

MDF-86V728E

MDF-86V838E

MDF-86V408E

MDF-86V340E

MDF-86V188E

MDF-86V58

58 -838 liters
2,400-60,000 samples

EXTERNAL STRUCTURE

File box
Easy to place log files

Ergonomic handle
Easy open
Lock&keys

Quick-release filter
Easy open & Quick clean

7" touch screen
Vaccum release port
USB data logger

Heavy duty castors
Universal wheels

PM ULT PM ULT freezers are made of color spayed anti-bacterial and anti-corrosion cold-rolled steel. Various configurations and humanized designs bring great convenience to users.

Access port


Remote interface
RS232/485

USB interface
Power on/off

FEATURERS

7" LCD touch screen
Easy operation and observation

Smart Operating System
Authorities, alarms, curve & events

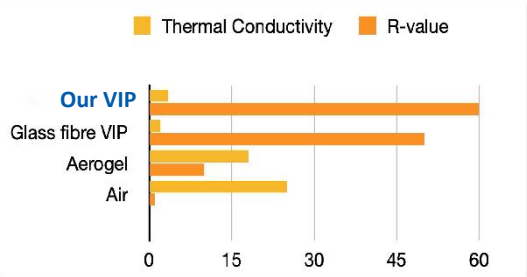

Advanced insulation tech
Fumed Silicon V.I.P & PURF

Multi-layer sealing design
Foamed inner doors with silicon gaskets

Eco-friendly
CFC-free refrigerant & insulation

Energy saving
More efficient cooling system



INTERNAL STRUCTURE




Adjustable shelves



PTC sensor



Pressure-regulated inner door handle

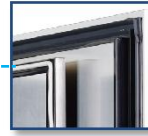

Multi-layer silicone sealing gaskets
Low temperature resistance



SUS 304 inner door filled with PURF
Detachable design

ALL
SUS 304
inside



SPECIFICATIONS

Basic		Refrigeration System		Alarms	
Cooling method	Direct Cooling	Refrigerant	R4001(mixed gas)	High/low temperature	Y
Defrost	Manual	Temp. controller	Microprocessor	Power failure	Y
Temperature range	-40℃~-86℃	System	Auto-cascade	Abnormal voltage	Y
Ambient temperature	10~32℃	Sensor	PTC 100	Sensor error	Y
Display	LED digital	Compressor brand/QTY	SECOP/2	Door ajar	Y
Access port QTY/Dia.	1/25mm	Capacity(L/cu.ft)	408/14.40	Filter screen check	Y
USB Intereface	Yes/1	Insulation	PURF+VIP	Thermostat failure	Y
Remote alarm Access	RS232/485	Inner door material/QTY	Foamed SUS304/2	Low battery	Y
Outside material	Color sprayed steel	Racks	8(5X4)+8(4X4)	Condenser clean	Y
Inside material	SUS304	Boxes (2 Inches)	288	High ambient temperature	Y
Castor	Universal	Samples	28,800	Power failure backup (alarm)	72h

Electrical Data		Dimensions		Accessories	
Voltage frequency(V/Hz)	220V 50Hz/60Hz	Package size(W*D*H)(mm)	954*955*2094	Chart recorder	Optional
Power(W)	890W	Exterior size(W*D*H)(mm)	872*1028*1945.5	Data logger	Optional
Current(A)	5.34A	Interior size(W*D*H)(mm)	606*575*1180	CO2 back up system	Optional
Power consumption	12.82KWh/24h	Shipping CBM	1.92	LN2 back up system	Optional
Noise(db)	53	N.W/G.W(kg)	268/293	Gloves	Optional

ACCESSORIES



LN₂/CO₂ backup system



Data logger



Chart recorder



Racks & Boxes



Gloves